

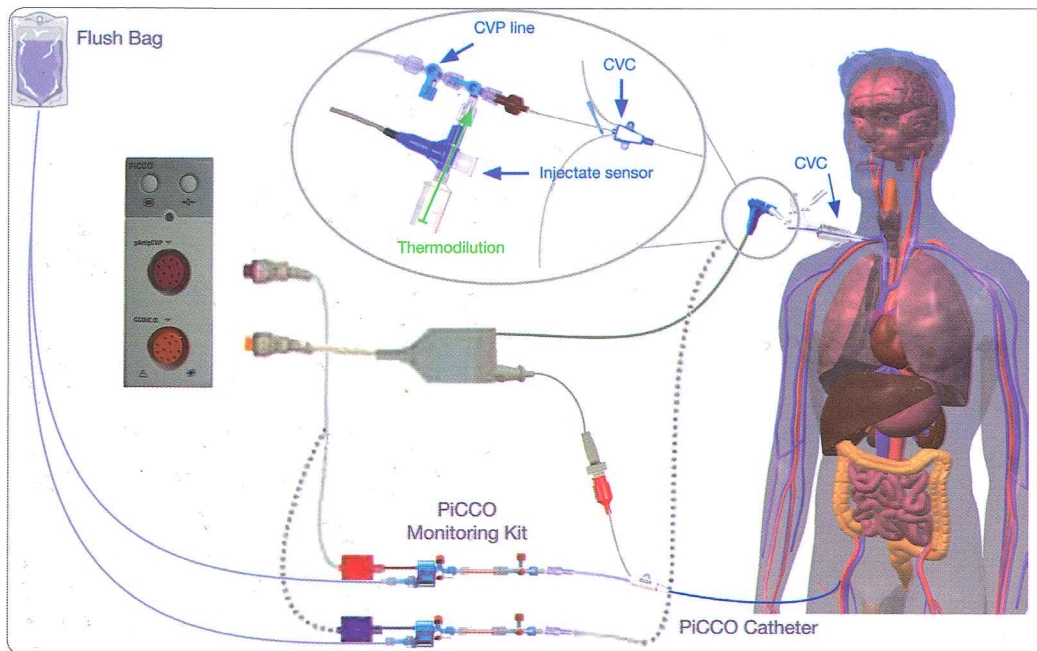


# Mindray PiCCO Setup Guide

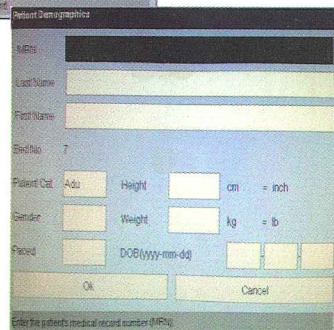

## 1. Setup

- Select **<Patient Setup>** from the tab on the main screen
- To admit a new patient select **<Admit Patient>** and input details
- For existing patients, select **<Patient Demographics>** and ensure that the following are entered;
  - Last Name • Patient Category • Height
  - Gender • Actual Weight • DOB
- Complete by selecting **<OK>**

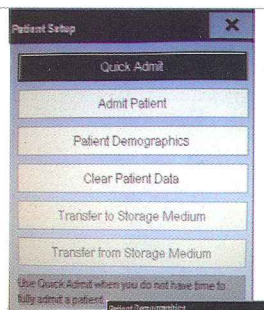


## 2. Patient Admission

# Mindray PiCCO Setup Guide

## 3. Screen Setup

- To activate PiCCO cardiac monitoring, select the waveform icon on the monitor
- On bottom left of screen there are 4 boxes. Touch one that is labelled OFF and select CCO
- To display a PiCCO monitoring screen, select the **<Choose Screen>** tab at the top of the display and select PiCCO
- Exit screen by touching **X** on top right

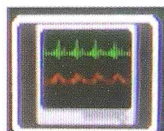

## 4. Zero Adjustment

- Perform zero adjustment of the PiCCO arterial line and central line as usual
- Open three way stop-cocks of the arterial and venous pressure transducers to the atmosphere
- For zero adjustment press  $\leftrightarrow 0 \leftrightarrow$  on the module or monitor
- Repeat zero adjustment prior to each set of thermodilutions or when necessary





- To setup for a set of thermodilutions, select **<Start PiCCO>** and then select **<Setup>**
- Ensure that patient's demographics have been entered correctly
- Set injectate volume to desired amount (15ml for patients <100kg & 20ml for those >100kg)
- Specify the PiCCO catheter position (femoral/ brachial/ axillary)
- Set CVP Measure to auto when measuring CVP continuously


*NB: C.O. Measure set to Auto will automatically perform a series of thermodilutions one after the other*



CCO Setup				
Height	175.0	cm	Inf. Volume	15ml
Weight	110.0	kg	Cat Type	PV2013L07
Patient Cat	Adu		Cat Position	Axillary/Brachial a
Gender	Male		C.O. Measure	Auto
PBW	70.6	kg	Exit PiCCO Screen	
BSA	1.900	m <sup>2</sup>	PiCCO Guide >>	
PBSA	1.900	m <sup>2</sup>	Select Parameter >>	
pCVP Measure	Auto		Hemo Para >>	
pCVP	9.5	cmH <sub>2</sub> O	Alarm Setup>>	
Enter the patient's height.				


- To ensure that there are no connection errors, check that the blood temperature (TB) is displayed and accurate (+/- X°C), and that the injectate temperature (TI) is ambient
- To begin thermodilution press **<Start>**
- Follow the instruction given on the screen, and when prompted; Inject the cold bolus of saline as quickly and smoothly as possible (< 7 seconds)
- Repeat this procedure until three agreeable results are obtained [see next section]

*NB: To increase accuracy use colder saline, inject faster or more smoothly. In patients with elevated lung water or very weak bloodflow; use 20ml instead of 15ml*

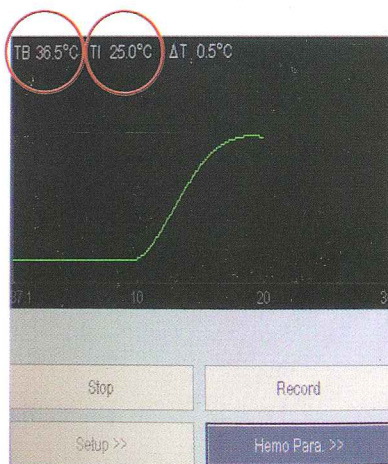

## 7. TD Quality


- ① A thermodilution curve must have a change in temperature ( $\Delta T$ ) of 0.20°C to be valid
  - ② The temperature of the injectate should be less than 10°C
  - ③ Three measurements in close agreement (+/- 10%) should be taken to ensure accuracy.
  - ④ Inaccurate measurements should be discarded by un-ticking the box on the left hand side
- Exit screen by touching **X** on top right of window

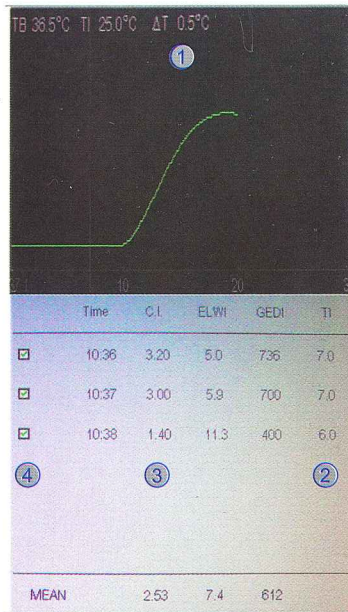

NB: Changes in temperature ( $\Delta T$ ) < 0.2°C may be due to temperature of saline, speed of injection or a patient condition (e.g. very weak blood flow, valve regurgitation, elevated lung water)

## 8. Visual Display

- To specify parameters to display, select **<PiCCO Start>-<Setup>-<Select Parameters>**. Here you can choose 3 parameters to monitor
- Spider Tab displays continuous PiCCO parameters and can be configured to user preference
- **<Haemodynamic Parameters>** tab displays all parameters calculated.
- **<Record>** will send to programmed printer.

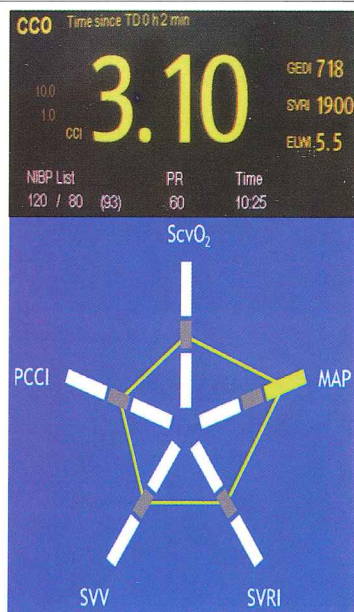



# Haemodynamic Decision Model

This decision model is not obligatory. It cannot replace the individual therapeutic decisions of the treating physician.

CI (l/min/m<sup>2</sup>)

Measured Values

GEDI (ml/m<sup>2</sup>)  
or ITBI (ml/m<sup>2</sup>)

ELWI (ml/kg)

Therapy Options

Targeted Values

1. GEDI (ml/m <sup>2</sup> ) or ITBI (ml/m <sup>2</sup> )	> 700 > 850	700-800 850-1000	> 700 > 850	700-800 850-1000
2. Optimise SVV (%)*	< 10	< 10	< 10	< 10
GEF (%)	> 25	> 30	> 25	> 30
or CFI (1/min)	> 4.5	> 5.5	> 4.5	> 5.5
ELWI (ml/kg) (slow response)		≤ 10		≤ 10

V+ = volume loading

V- = volume reduction

Cat = catecholamine / cardiovascular agents

\*SVV is only applicable in fully ventilated patients without cardiac arrhythmia

**PULSION**  
Medical Systems

PULSION Medical Systems SE • Hans-Riedl-Str. 17 • D-85622 Feldkirchen, Germany  
Tel. +49-(0)89-45 99 14-0 • Fax +49-(0)89-45 99 14-18  
info@pulsion.com • www.pulsion.com

CE 0124

## Haemodynamic Normal Values

### Central Venous Oxygenation - Oxygenation Balance

(Oxygen load of the venous blood after passing through the organs)

ScvO<sub>2</sub> \*\* 70-80 %

O<sub>2</sub> Consumption (Consumption of O<sub>2</sub> by organs)

VO<sub>2</sub>I 125-175 ml/min/m<sup>2</sup>

O<sub>2</sub> Delivery (Delivery of O<sub>2</sub> via blood to organs)

DO<sub>2</sub>I 400-650 ml/min/m<sup>2</sup>

Haemoglobin (Oxygen transporter in blood)

Hb \*\*\* 8.7-11.2 mmol/l (Male)  
7.5-9.9 mmol/l (Female)

Arterial / capillary oxygen saturation (Oxygen load of arterial blood)

SaO<sub>2</sub> / SpO<sub>2</sub> 96-100 %

Flow

Cardiac Index

CI

3-5 l/min/m<sup>2</sup>

Pulse Contour Cardiac Index (Cardiac Index related to body surface)

PCCI

3-5 l/min/m<sup>2</sup>

Chronotropy

Heart Rate

HR

60-80 bpm

Stroke Volume

Preload

Stroke Volume Index (Output per heart beat)

SVI

40-60 ml/m<sup>2</sup>

Global Enddiastolic Volume Index (Volume of blood in the heart)

GEDI

680-800 ml/m<sup>2</sup>

Intrathoracic Blood Volume Index (Volume of blood in heart and lungs)

ITBI

850-1000 ml/m<sup>2</sup>

Stroke Volume Variation (Dynamic fluid responsiveness)

SVV \*

0-10 %

Pulse Pressure Variation (Dynamic fluid responsiveness)

PPV \*

0-13 %

Afterload

Systemic Vascular Resistance Index (Resistance of vascular system)

SVRI

1700-2400 dyn\*sec/cm<sup>5</sup>\*m<sup>2</sup>

Mean Arterial Pressure

MAP

70-90 mmHg

Contractility

Global Ejection Fraction (Ratio of stroke volume and preload)

GEF

25-35 %

Left Ventricular Contractility (Increase of arterial pressure over time)

dPmax

Trend information

Cardiac Function Index (Ratio of CI and preload)

CFI

4.5-6.5 1/min

Cardiac Power Index (Global cardiac performance)

CPI

0.5-0.7 W/m<sup>2</sup>

Lung

Extravascular Lung Water Index (Lung oedema)

ELWI

3-7 ml/kg

Pulmonary Vascular Permeability Index (Permeability of lung tissue)

PVPI

1.0-3.0

Liver

Plasma Disappearance Rate ICG (Performance of the liver)

PDR

16-25 %/min

Retention rate of ICG after 15 minutes (Performance of the liver)

R15

0-10 %

Absolute values (non-indexed values) are only usable in trend screens and have no normal range. \* SVV/PPV are only applicable in fully ventilated patients without cardiac arrhythmias.

\*\* A high-normal / high ScvO<sub>2</sub> can be a sign of insufficient O<sub>2</sub> utilisation \*\*\* 14-18 g/dl (Male); 12-16 g/dl (Female)