



Chest Drainage System

10 Easy Steps To Replace A Used Bottle

Rocketmedical

REPLACEMENT BOTTLE SETUP



- 1 Rotate red sealing cap to open.
- 2 Fill to prime level with sterile water or normal saline, approx 500ml.




- 5 Insert tubing set into new bottle and twist Clockwise until 'click' is heard.
- 6 Unclamp catheter.




If suction is used

- 7 Remove green vent cap.
- 8 Connect suction tubing to vent/suction port.

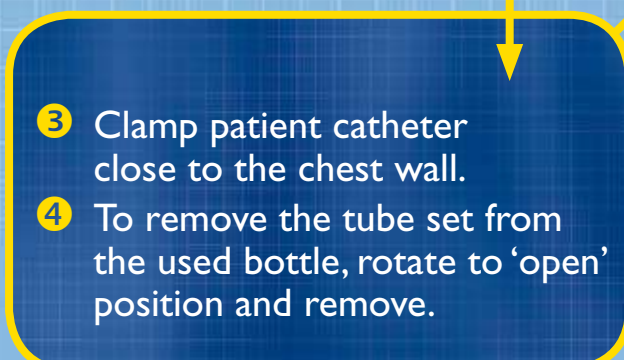

- 3 Clamp patient catheter close to the chest wall.
- 4 To remove the tube set from the used bottle, rotate to 'open' position and remove.

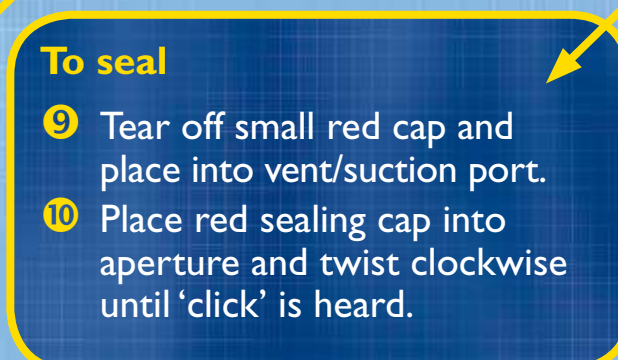

- ### To seal
- 9 Tear off small red cap and place into vent/suction port.
 - 10 Place red sealing cap into aperture and twist clockwise until 'click' is heard.



Used bottle is now ready for disposal. According to hospital procedure.

RE-ORDER CODES			
Single Chamber			
Blue Bottle	R54.500	Floor Stand	R54.503
Tubing Set		Double bottle	
-with connector	R54.502	Floor Stand	R54.505
-without connector	R54.539	Dual Collection Bottle	
		Green Bottle	R54.509

USED BOTTLE DISPOSAL



UK Customer Services

Rocket Medical plc Sedling Road Washington Tyne & Wear NE38 9BZ England
Telephone: +44 (0) 191 419 6988 Fax: +44 (0) 191 419 6989
Email: homesales@rocketmedical.com

Head Office & Export Customer Services

Rocket Medical plc Imperial Way Watford Herts WD24 4XX England
Telephone: +44 (0) 1923 651404 Fax: +44 (0) 1923 240334
Email: exportsales@rocketmedical.com