


## Intra-Abdominal Pressure Measurement

### Equipment Needed:

- Pressure Module, Cable and Transducer
- 500ml bag 0.9% Saline and IV giving set
- 60ml luer lock syringe
- 3 way taps x2
- Green 21g needle

### Assemble Equipment as Below:

- Zero transducer at mid axillary line the same as vascular pressure lines

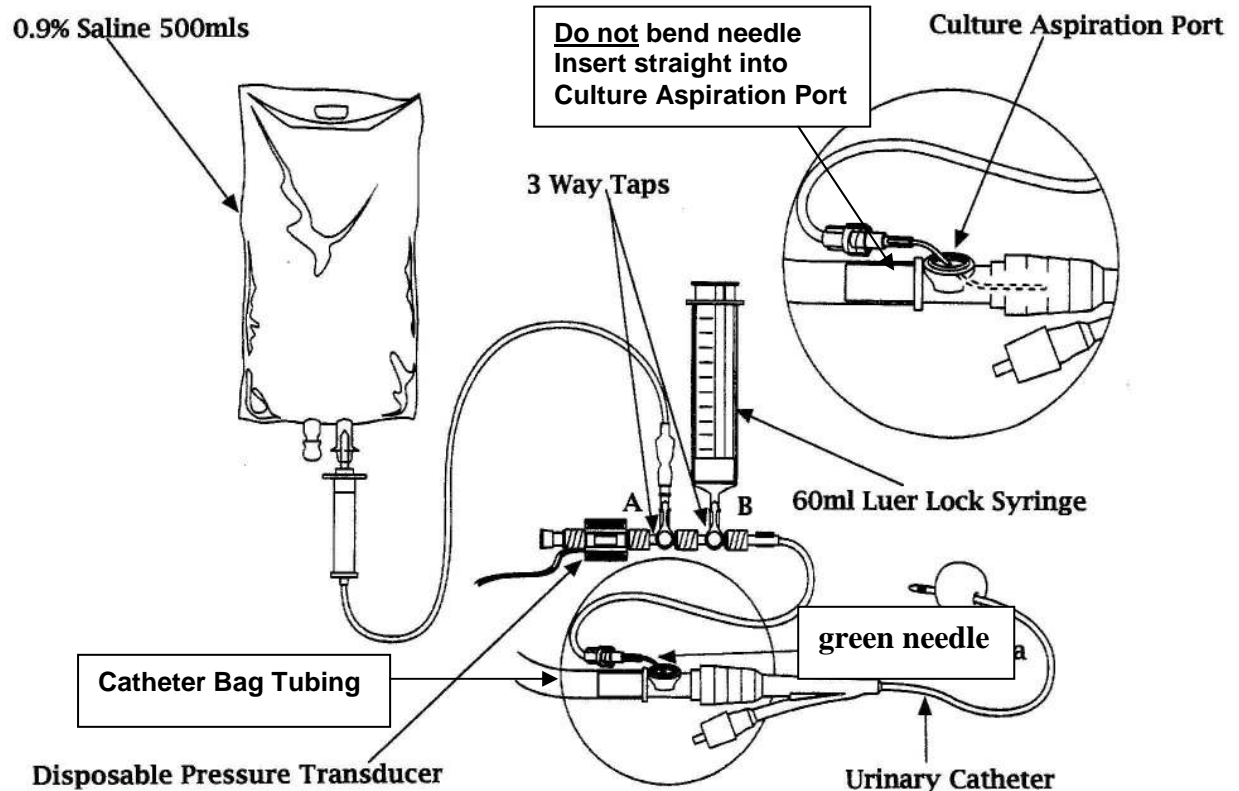

### MEASUREMENT:

1. Clamp catheter bag tubing
2. Tap B OFF to urine catheter, Tap A OFF to pressure transducer
3. 50ml Saline aspirated into syringe from bag and turn Tap A OFF to bag
4. Tap B OFF to pressure transducer and inject 50ml saline into bladder
5. Tap B OFF to syringe
6. Catheter bag tubing clamp- release momentarily to allow all air to flushout and reclamp
7. Read intra-abdominal pressure from monitor (end expiration) and record on chart
8. Catheter bag tubing clamp- release
9. Deduct 50ml (saline volume instilled) from next hour's urine output
10. Disconnect green needle and discard into sharps bin

<b>Title:</b> Intra-abdominal pressure measurement	
<b>ID:</b>	<b>Authors</b> B Cook
<b>Status Draft/Final:</b> Final	<b>Approved by:</b> QIT editorial group
	<b>Written:</b> June 2004
<b>Reviewed:</b> March 2009, Nov 2012, 2019	<b>Review Date:</b> June 2021